

Connecting your Nathos Nova M to the myPhonak app

Download the myPhonak app


- 1 Make sure Bluetooth® is on
- 2 Install the myPhonak App from Google Play™ or Apple App Store®

Make sure your smartphone is connected to the internet

- 3 In order to use the app, you must tap on I agree to accept the Privacy Notice to continue.
- 4 Product improvement: You can optionally share usage data, which enables us to learn and improve our products.


Enter pairing mode

Open and close the battery door, you will then have 3 minutes to complete the pairing.

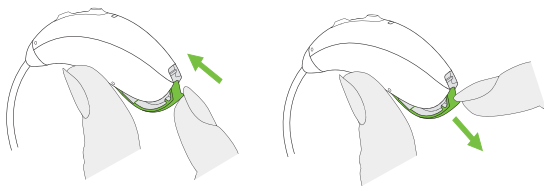

Pairing the hearing aids to the myPhonak app


1

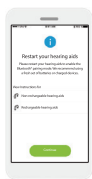

Pairing instructions


Tap Continue to initiate the search process. Tap the instructions for non-rechargeable hearing aids to review the instructions for your device.

2

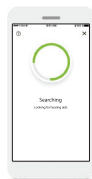


Searching

The app is searching for compatible hearing aids and will display them once they are detected. This may take a while.

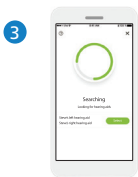

3

Selecting
Tap Select when your hearing aids appear in the list.

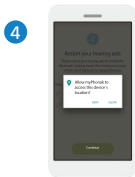

4

Location
On Android devices, you must enable location services when pairing Bluetooth devices for the first time. After the initial setup you can disable the location services again.



5

Pairing the hearing aids
The app will connect to each hearing aid separately.



6

Confirm for all hearing aids
Please confirm by tapping Pair in the pop-up for every device separately.



7

Pairing complete
Both hearing aids are now paired. The app will automatically proceed to the next step.



8

Setup complete
You are now ready to use the myPhonak app. Tap OK to access the main screen.